



HTL PIPETTE TIPS



***perfect seal between tip
and pipette shaft***

***accurate graduations for visual
confirmation of sample dispensed***

***optimized bevel for the most
accurate dispensing***


***DNase, RNase and endotoxin-free
for highest sample integrity***

***filter in a tip prevents aerosol
liquid contamination***

Ordering information

Bulk: 1000 tips/pack, 20 packs/case
Racked: 96 tips/rack, 10 racks/unit, 5 units/case

Stack Rack: 5 x 96 tips/stack, 10 stacks/case
Refill System: 2 stacks of 5 x 96 tips/unit, 5 units/case

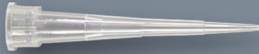Volume	Tip type	Packaging	Cat. no.
0.5 - 10 μ l	Microvolume	Bulk Racked Racked, Pre-Sterilized Stack Rack Refill System	30011 30012 30013 30014 30015
 Fits: LM/D/DV2, LM/D/DV10, DV8-10, DV12-10			
1 - 200 μ l	Bevelled	Bulk Racked Racked, Pre-Sterilized Stack Rack	30201 30202 30203 30204
 Fits: LM/D/DV20, LM/D/DV50, LM/D/DV100, LM/D/DV200, DV8-50, DV8-200, DV12-50, DV12-200			
1 - 200 μ l Graduation marks at 10, 50 and 100 μ l	Bevelled, Reference	Bulk Racked Racked, Pre-Sterilized Stack Rack Refill System	30221 30222 30223 30224 30225
 Fits: LM/D/DV20, LM/D/DV50, LM/D/DV100, LM/D/DV200, DV8-50, DV8-200, DV12-50, DV12-200			
1 - 300 μ l	Fine-Point	Bulk* Racked Racked, Pre-Sterilized	30301 30302 30303
 Fits: LM/D/DV250, DV8-300, DV12-300		Bulk: 1000 tips/pack, 10 packs/case	
100 - 1000 μ l	Bevelled	Bulk* Racked* Racked*, Pre-Sterilized	31001 31002 31003
 Fits: LM/D/DV1000		Bulk: 1000 tips/pack, 5 packs/case Racked: 100 tips/rack, 10 racks/unit, 5 units/case	
100 - 1200 μ l	Bevelled	Bulk* Racked* Racked*, Pre-Sterilized	31201 31202 31203
 Fits: LM/D/DV1000		Bulk: 1000 tips/pack, 5 packs/case Racked: 96 tips/rack, 12 racks/unit, 3 units/case	
5 ml	Macrovolume	Bulk*	35001
 Fits: LM/D/DV5000		Bulk: 250 tips/pack, 10 packs/case	
10 ml	Macrovolume	Bulk*	36001
 Fits: LM/D/DV10000		Bulk: 200 tips/pack, 10 packs/case	

Filter tips

Racked: 96 tips/rack, 10 racks/unit, 5 units/case

Volume	Tip type	Packaging	Cat. no.
0.5 - 10 µl	Microvolume	Racked, Pre-Sterilized	30016
			
Fits: LM/D/DV2, LM/D/DV10, DV8-10, DV12-10			
0.5 - 20 µl	Universal	Racked, Pre-Sterilized	30026
			
Fits: LM/D/DV20			
0.5 - 50 µl	Universal	Racked, Pre-Sterilized	30056
			
Fits: LM/D/DV50, DV8-50, DV12-50			
1 - 100 µl	Universal	Racked, Pre-Sterilized	30106
			
Fits: LM/D/DV100			
1 - 200 µl	Universal	Racked, Pre-Sterilized	30206
			
Fits: D/DV200, DV8-200, DV12-200			
1 - 300 µl	Universal	Racked, Pre-Sterilized	30306
			
Fits: LM/D/DV250, DV8-300, DV12-300			
100 - 1000 µl	Universal	Racked*, Pre-Sterilized	31006
			
Fits: LM/D/DV1000		Racked: 100 tips/rack, 10 racks/unit, 5 units/case	
100 - 1100 µl	Bevelled	Racked*, Pre-Sterilized	31206
			
Fits: LM/D/DV1000		Racked: 96 tips/rack, 12 racks/unit, 3 units/case	



HTL Pipette Tips are designed to fit a wide variety of single and multichannel pipettes. Flexible walls and the diameter of tips collars ensure a secure fit with less force required to load and eject. Aerosol barrier tips prevent aerosol and liquid contamination. Technologically advanced filter material using in pipette tips allows free passage of air while preventing passage of aerosols. Filters contain no leachable additives. All pre-sterilized tips are exposed to e-beam radiation for complete security without degrading the plastic or performance.

Tip packaging options




BULK PACKAGING

Pipette tips are available in bulks - resealable plastic bags - keeping them safe from contamination. One bulk contains 200, 250 or 1000 pieces of tips depends on the tip volume.

HINGED RACKS

Tip racks in durable polypropylene provide good stability on the lab bench. They feature a hinged lid with a snap fastener. One rack contains 96 or 100 pieces of tips depends on the tip volume.

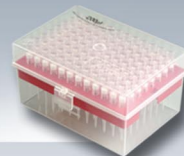

STACK RACK

Inter-locking stack racks offer stability for the tips while saving valuable lab bench space. One stack rack contains 5 trays with 96 tips. Stacks are available only for 10 μ l and 200 μ l tips.

REFILL SYSTEM

Space saving refill stacks for easy reload of an empty rack. One refill contains 10 trays with 96 tips each. Refills are available only for 10 μ l tips and 200 μ l tips with graduation marks.



Autoclaving guidelines

To ensure successful steam sterilization of your tips please adhere to the following guideline: Autoclave all tips at 121°C for 15 minutes, 100 kPa. Do not use dry cycle as this may cause deformation or weakening of the plastic. Allow the temperature within the autoclave to return to at least 80°C before removing product.

Distributed by:

